

TEST REPORT

(Electronic version)



No: 200092748

VERIFICATION WEBSITE: www.gttc.net.cn

VERIFICATION CODE: KABB-4733-24

ISSUE DATE: 2020-04-23

二维码下载报告



APPLICANT: GUANGDONG FOR YOU ESSENTIALS MANUFACTURING CO., LTD
ADDRESS: NO. 6, SHANHUA ROAD, YANGHE TOWN, GAOMING DISTRICT, FOSHAN CITY, GUANGDONG, CHINA

INFORMATION CONFIRMED BY APPLICANT:

KN95 PROTECTIVE MASK

QUANTITY: THIRTY PIECES

BRAND: V-SHINE

COLOUR: WHITE

MANUFACTURE'S NAME: GUANGDONG FOR YOU ESSENTIALS MANUFACTURING CO., LTD

DATE RECEIVED/DATE TEST STARTED: 2020-04-22

CONCLUSION:

FILTRATION EFFICIENCY TO NaCl PARTICULATE MATTER	M
INSPIRATORY RESISTANCE	M
EXPIRATORY RESISTANCE	M
FLAMMABILITY	M
HEAD BAND	M

NOTE: "M" -MEET THE STANDARD'S REQUIREMENT "F" -FAIL TO MEET THE STANDARD'S REQUIREMENT
"----" -NO COMMENT

REMARK:

THIS REPORT IS THE ENGLISH TRANSLATION VERSION OF THE REPORT 200059110.
ALL THE TESTED ITEMS ARE TESTED UNDER THE STANDARD CONDITION (EXCEPT FOR INDICATION).
COPIES OF THE REPORT ARE VALID ONLY RE-STAMPED.
THE EXPERIMENT WAS CARRIED OUT AT No. 1, ZHUJIANG ROAD, PANYU DISTRICT, GUANGZHOU, GUANGDONG, P. R. CHINA.

APPROVED BY:

ZiShan Guo SENIOR ENGINEER

郭子山




PAGE 1 OF 4

TEST REPORT

(Electronic version)

No:200092748

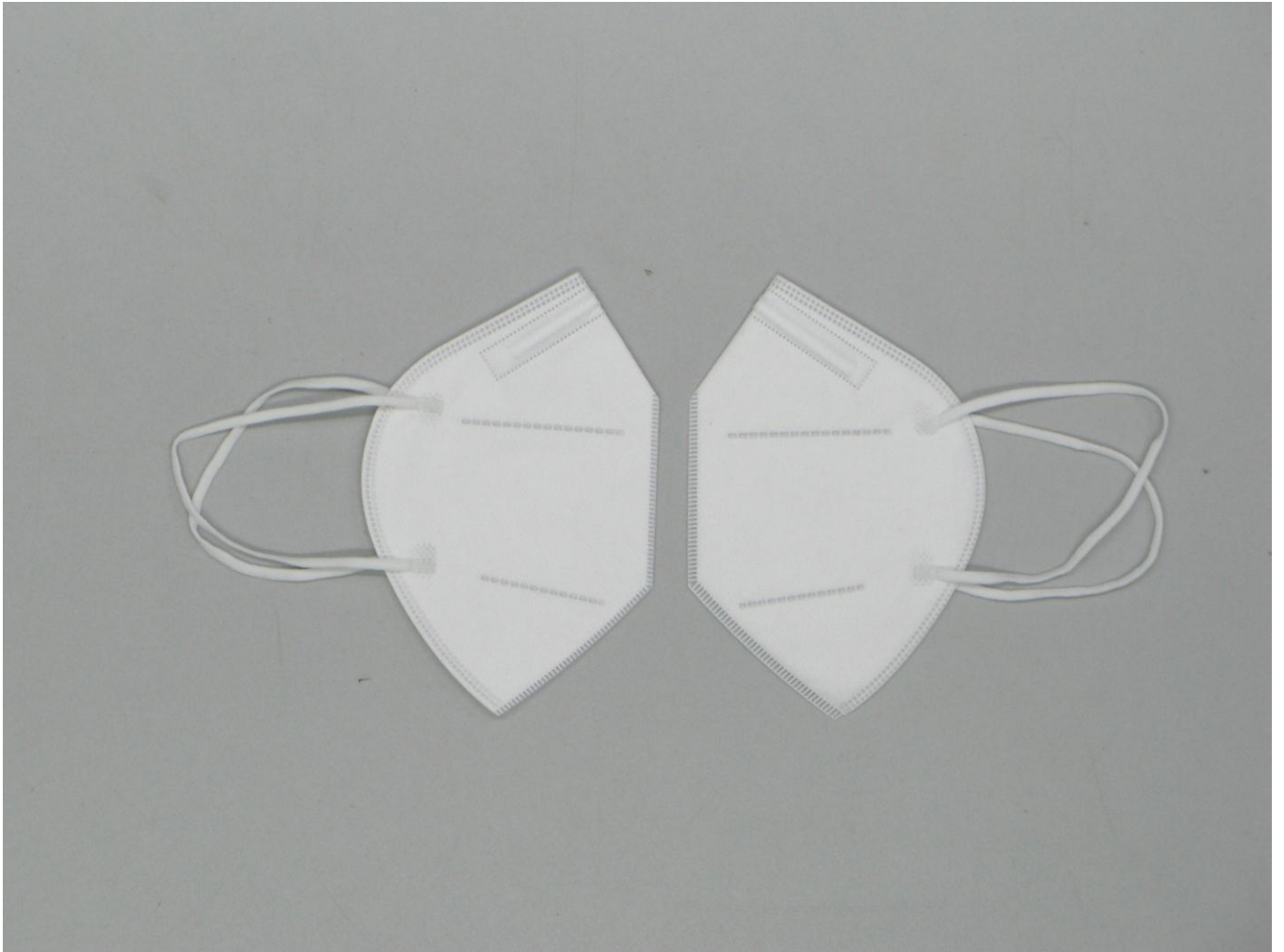

TEST REPORT

(Electronic version)

No:200092748

FILTRATION EFFICIENCY TO NaCl PARTICULATE MATTER

(GB 2626-2006 6.3, AIR FLOW:85L/min, AEROSOL:NaCl, AEROSOL CONCENTRATION:15mg/m³, TEMP:23.1℃, RH:36.2%)

FILTRATION EFFICIENCY(%) :

UNTREATED SAMPLE

1# 99.666
2# 99.689
3# 99.811
4# 99.584
5# 99.612
6# 99.720
7# 99.675
8# 99.718
9# 99.650
10# 99.628

CONDITIONING TREATED

1# 99.736
2# 99.578
3# 99.624
4# 99.611
5# 99.780

REQUIREMENT

FILTRATION EFFICIENCY(%) :

≥95.0
(KN95)
(GB 2626-2006)

INSPIRATORY RESISTANCE (Pa)

(GB 2626-2006 6.5, HEAD SIZE: MEDIUM)

UNTREATED SAMPLE:

1# 94.5
2# 96.3

PRETREATMENT SAMPLE:

1# 87.5
2# 85.8

REQUIREMENT

≤350
(GB 2626-2006)

EXPIRATORY RESISTANCE (Pa)

(GB 2626-2006 6.6, HEAD SIZE: MEDIUM)

UNTREATED SAMPLE:

1# 61.0
2# 63.2

PRETREATMENT SAMPLE:

1# 59.9
2# 58.6

REQUIREMENT

≤250
(GB 2626-2006)

FLAMMABILITY (s)

(GB 2626-2006 6.15)

AFTERFLAME TIME

UNTREATED SAMPLE

1# 0.0
2# 0.0

CONDITIONING TREATED

3# 0.0
4# 0.0

REQUIREMENT

AFTERFLAME TIME
≤5
(GB 2626-2006)



PAGE 3 OF 4

TEST REPORT

(Electronic version)

No:200092748

HEAD BAND

(GB 2626-2006 6.11)

UNTREATED SAMPLE:

1# PASS

CONDITIONING TREATED:

1# PASS

REQUIREMENT

EACH HEAD BAND, BUCKLE AND OTHER ADJUST
PARTS OF MASK SHOULD NOT SLIP OR BREAK
UNDER 10N TENSION FOR 10S.
(GB 2626-2006)



—End of Report—

PAGE 4 OF 4