

Test Report **SL52025233245001TX** **Date: April 22, 2020** **Page 1 of 3**
GUANGDONG FOR YOU ESSENTIALS MANUFACTURING CO.,LTD
NO.6, SHANHUA ROAD, YANGHE TOWN, GAOMING DISTRICT, FOSHAN, CHINA

The following sample(s) was/were submitted and identified on behalf of the client as:

Sample Description : (A)Mask


Style No. : 17.5*9.5CM,3 layers
Composition : (A)Melt blown MB / non-woven PP
Sample Color : (A)blue

Proposed Care Instruction : /

Test Performed : Selected test(s) as requested by applicant

Sample Receiving Date : Mar 23,2020
Testing Period : Mar 31,2020 - Apr 22,2020
Test Result(s) : For further details, please refer to the following page(s).

Signed for and on behalf of
SGS-CSTC Standards Technical Services (Shanghai) Co., Ltd Testing Center

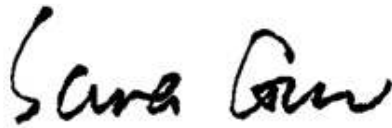

Sara Guo (Account Executive)

Unless otherwise agreed in writing, this document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law. Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.

Attention: To check the authenticity of testing /inspection report & certificate, please contact us at telephone: (86-755) 8307 1443, or email: CN.Doccheck@sgs.com



SGS-CSTC Standards Technical Services (Shanghai) Co., Ltd.
Testing Center

3rd Building, No.889, Yishan Road, Xuhui District Shanghai, China 200233
中国·上海·徐汇区宜山路889号3号楼 邮编: 200233

t (86-21) 61402666 f (86-21) 64958763
t (86-21) 61402666 f (86-21) 64958763

www.sgsgroup.com.cn
e sgs.china@sgs.com

Test Result

Medical Face Masks-Requirements and Test Methods

(EN 14683:2019)

Clause 5.2.2 Bacterial filtration efficiency (BFE)@**

| | 1# | 2# | 3# | 4# | 5# |
|----------|------|------|------|------|------|
| (BFE), % | 99.9 | 99.7 | 99.7 | 99.8 | 99.8 |

Remark: Performance Requirement: Type I>95%, Type II>98%, Type IIR >98%

** : The test was carried out by external laboratory assessed as competent

@ : This test method is not in CNAS accredited scope

Clause 5.2.3 Breathability

(EN 14683 :2019 Annex C, Flow rate 8 l/min)

Sample: A

| | 1# | 2# | 3# | 4# | 5# |
|---|------|------|------|------|------|
| Differential pressure ΔP (Pa/cm ²) | 42.3 | 46.5 | 58.3 | 47.1 | 47.0 |

Remark: Performance Requirement: Type I<40 Pa/cm², Type II<40 Pa/cm², Type IIR<60 Pa/cm²

Clause 5.2.4 Splash Resistance

(ISO 22609 :2004, Pressure 16.0 kPa)

Sample: A

Penetration on inside surface

| | | | | | | | |
|------|------|------|------|------|------|------|------|
| 1# | 2# | 3# | 4# | 5# | 6# | 7# | 8# |
| Pass | Pass | Pass | Pass | Pass | Pass | Pass | Pass |
| 9# | 10# | 11# | 12# | 13# | 14# | 15# | 16# |
| Pass | Pass | Pass | Pass | Pass | Pass | Pass | Pass |
| 17# | 18# | 19# | 20# | 21# | 22# | 23# | 24# |
| Pass | Pass | Pass | Pass | Pass | Pass | Pass | Pass |
| 25# | 26# | 27# | 28# | 29# | 30# | 31# | 32# |
| Pass | Pass | Pass | Pass | Pass | Pass | Pass | Pass |

Number of Pass: 32
Overall result: Acceptable

Remark:

- 1) Performance Requirement Type I: N/A, Type II: N/A, Type IIR: ≥16.0kPa
- 2) Distance of the medical face mask target area surface to the tip of cannula is 300±10mm.
- 3) Condition and Test temperature (21±5)°C, relative humidity (85±10)%
- 4) An acceptable quality limit of 4.0% is met for a single sampling plan when 29 or more of the 32 tested specimens show pass results

Unless otherwise agreed in writing, this document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law. Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.

Attention: To check the authenticity of testing /inspection report & certificate, please contact us at telephone: (86-755) 8307 1443, or email: CN.Doccheck@sgs.com



Clause 5.2.5 Microbial Cleanliness
(EN 14683: 2019 Annex D)

| | 1# | 2# | 3# | 4# | 5# |
|-------|----|----|----|----|----|
| CFU/g | 1 | 1 | 1 | <1 | <1 |

Remark: Performance Requirement: Type I ≤ 30 CFU/g, Type II ≤ 30 CFU/g, Type IIR ≤ 30 CFU/g

The statement of conformity in this test report is only based on measured values by the laboratory and does not take their uncertainties into consideration.

End of Report



Unless otherwise agreed in writing, this document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law. Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.

Attention: To check the authenticity of testing /inspection report & certificate, please contact us at telephone: (86-755) 8307 1443, or email: CN.Doccheck@sgs.com

SGS-CSTC Standards Technical Services (Shanghai) Co., Ltd. 3rd Building, No. 889, Yishan Road, Xuhui District Shanghai, China 200233 t (86-21) 61402666 f (86-21) 64958763 www.sgs.com.cn
 中国·上海·徐汇区宜山路889号3号楼 邮编: 200233 t (86-21) 61402666 f (86-21) 64958763 e sgs.china@sgs.com